


PRECISION SPINE  
**SURELOK™** MIS 3L  
PERCUTANEOUS SCREW SYSTEM

**SIMPLE & EFFICIENT**

The SureLOK™ MIS 3L is a low profile, percutaneous screw system. The system features an extended tab tulip with 20mm of controlled rod reduction. The screw has a triple lead, proximally tapered thread that greatly enhances delivery. This reduces fatigue and minimizes OR time without compromising screw purchase. In addition, the open tulip design, 30° of rod angulation and bulletted rod tip eases rod insertion. The system also offers multiple sleeve options to allow for aggressive tulip manipulation, rod reduction and compression and distraction techniques.

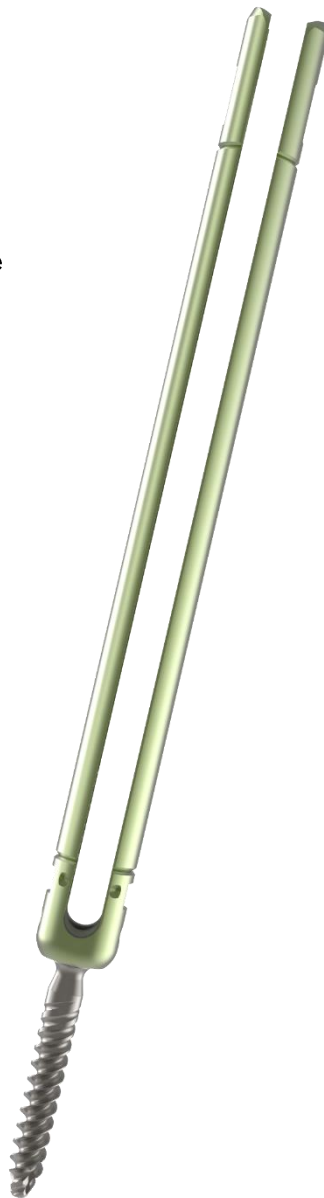


MULTIPLE SLEEVE OPTIONS



AGGRESSIVE SCREW TIP



TRIPLE LEAD THREAD

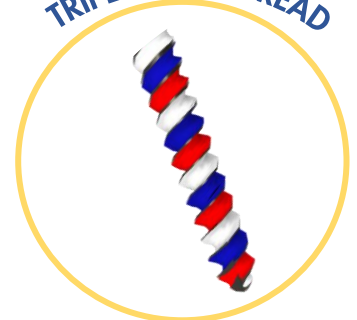

PRECISION SPINE®

*Discover the Difference*



# SURELOK™ MIS 3L PERCUTANEOUS PEDICLE SCREW

**150mm Integrated Extension Tab** features a low profile outer diameter (13.7mm) allowing for a smaller incision to minimize muscle disruption

**Open Tab Design** provides maximum visualization and eases rod placement

**Polyaxial Screws** offer 30° of angulation in all planes for intraoperative flexibility


**Integrated Reduction Threads** provide 20mm of controlled rod reduction

**Break-Off Tab** removal is completed in one easy step

**Triple Lead Thread** allows for efficient screw delivery

**Aggressive Self Tapping Screw Tip** optimizes screw engagement

**Curved & Straight Rod's** bullet tip and keyed hex end provide a rigid connection while the 15° rod inserter angle eases rod insertion and release

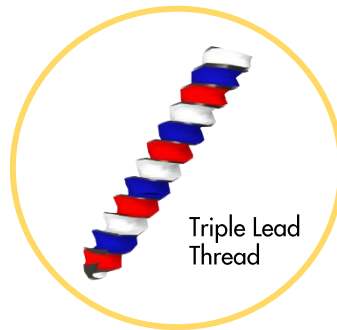

**Center and Side Load Rod Inserters** meet varying surgical techniques and provides a strong rod interface for secure insertion

Center Load Inserter



Side Load Inserter



**Streamlined Insertion Instrumentation** provides OR efficiency while minimizing soft tissue disruption:

- **Dilation/Awl/Tap** is completed in one step



- **Modular Rings & Sleeves** provide maximum rigidity



- **Internal Rod Reducer** allows for 50mm of controlled rod reduction



- **Compressor/Distractor** is powerful, precise and effective

## MIS 3L Rod Sizing – 5.5mm ø

Ti Curved	Ti Straight
35 - 80mm (5mm)	35-80mm (5mm)
90mm-110mm (10mm)	90-150mm (10mm)

## MIS 3L Sizing

5.5mm Diameter	6.5mm Diameter	7.5mm Diameter	8.5mm Diameter
35mm	35mm	35mm	35mm
40mm	40mm	40mm	40mm
45mm	45mm	45mm	45mm
50mm	50mm	50mm	50mm
55mm	55mm	55mm	55mm

### Precision Spine, Inc.

2050 Executive Drive, Pearl, MS 39208  
 Customer Service: 1.888.241.4773  
 Phone: 601.420.4244  
 Toll Free: 877.780.4370  
 Fax: 601.420.5501  
[www.precisionspineinc.com](http://www.precisionspineinc.com)

